

Storage chute

Model 4614, 4621, 6914 & 6921



Always available, always at glance

Made of lightweight, durable aluminum, our containers are designed for safe, organized, and hygienic storage of medical supplies such as band-aids, swabs, disposable syringes, tubes and other useful items. The powder-coated metal features high acid and temperature resistance. The smooth surfaces are easy and hygienic to clean.

The storage trays are swiveling and easy to fill. Thanks to the FIFO principle (First in - First Out), the material with the nearest material with the closest expiration date is always used. Filling takes place from above, removal from below. In this way, items that are filled first are also removed first. This means you avoid wasting material and minimize your disposal costs.

Details





Extraction





Open chute

Made in Germany

All of our models and the corresponding chutes are made in Germany.

 **Short supply chains**
 **CO₂-Savings**

 **No sea containers**
 **High quality**



Highlights

- For optimal storage of consumables
- Hygienic and easy to clean
- First In - First Out system (FIFO)
- Guaranteed order with immediate overview
- Individually adaptable to the needs of the respective workplace
- Flexible mounting and assembly options



Repainting possible

As a standard, the chutes are delivered uniformly in traffic white RAL 9016. Repainting according to RAL color code is possible.

Keyhole bore

Keyhole bore for mounting on the wall.

Vesa 75

VESA 75 standard for mounting on the wall.

Modular


The modules are available in different layouts and can be combined with each other.

Various sizes

The cases and internal dividers are available in different sizes.


FIFO

Extraction with closed chute: First In - First Out system.

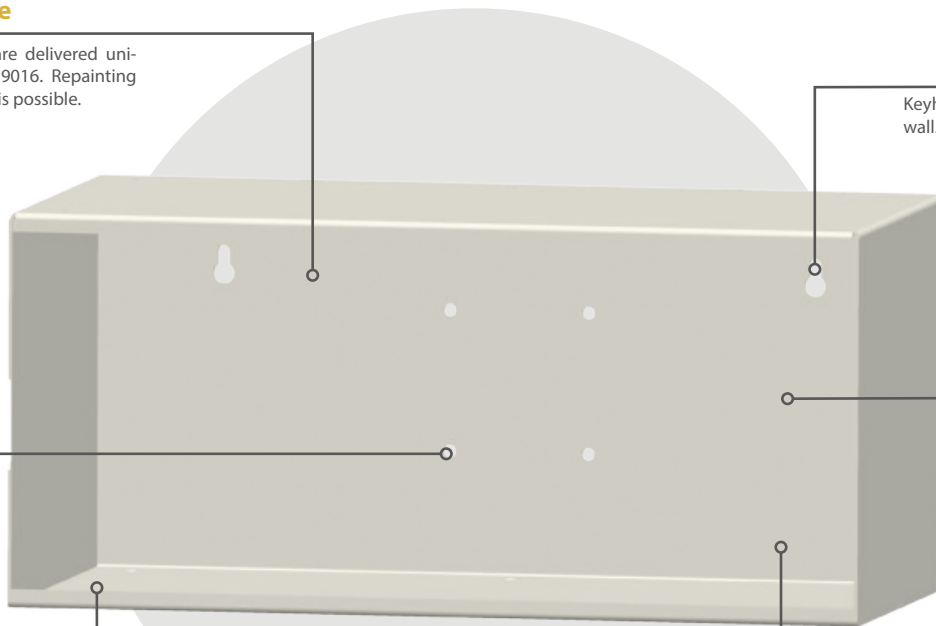

Flexible use due to modular design

Die Schüttenmodule sind in verschiedenen Größenausführungen mit unterschiedlichen Aufteilungen erhältlich und untereinander kombinierbar. Auch können Sie Ihre Wunschkombination vormontiert bei uns bestellen. Verschiedene Bohrungen in der Gehäuserückseite ermöglichen eine variable Befestigung. So können Sie die Schütten optimal an die Bedingungen und Ansprüche Ihres Arbeitsplatzes anpassen.

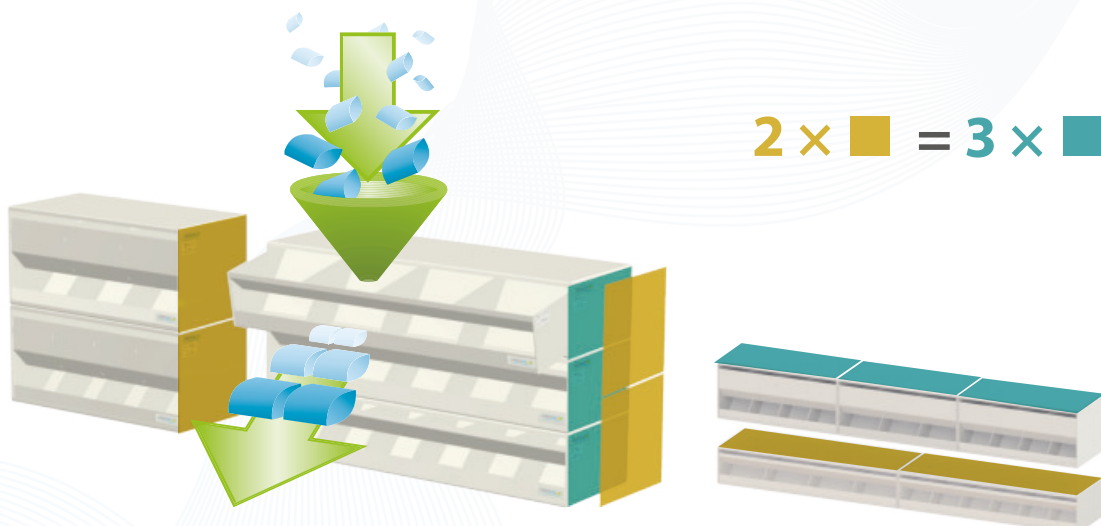

2 modules XX21 correspond to the height of 3 modules XX14 // 2 modules 46XX correspond to the width of 3 modules 69XX

Technical specifications

Size and weight:

Model 4614: 460 x 140 x 225 mm,
Weight incl. chute: ca. 4,5 kg
Model 4621: 460 x 210 x 225 mm,
Weight incl. chute: ca. 5,5 kg
Model 6914: 690 x 140 x 225 mm,
Weight incl. chute: ca. 5,5 kg

Model 6921: 690 x 210 x 225 mm,
Weight incl. chute: ca. 7 kg

Coloring:

Traffic white RAL 9016, repainting possible on request

Material:




Aluminum, powder coated

Sterility:

No

Order information





Model 4614

Art.-No.	Description	Image
311.1146140	Case module 4614	
311.1146143	Chute for module 4614, 3-fold division	
311.1146144	Chute for module 4614, 4-fold division	




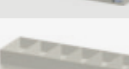
Model 4621

Art.-No.	Description	Image
311.1146210	Case module 4621	
311.1146213	Chute for module 4621, 3-fold division	
311.1146214	Chute for module 4621, 4-fold division	

Model 6914

Art.-No.	Description	Image
311.1169140	Case module 6914	
311.1169144	Chute for module 6914, 4-fold division	
311.1169145	Chute for module 6914, 5-fold division	
311.1169146	Chute for module 6914, 6-fold division	

Model 6921

Art.-No.	Description	Image
311.1169210	Case module 6921	
311.1169214	Chute for module 6921, 4-fold division	
311.1169215	Chute for module 6921, 5-fold division	
311.1169216	Chute for module 6921, 6-fold division	

Order and advice by:

HERZmed GmbH

Herkulan-Schwaiger-Gasse 16
82487 Oberammergau

Tel.: +49 (0) 911 / 14 88 959 - 0
Fax: +49 (0) 911 / 14 88 959 - 9

Mail: herzmed@herzmed.de
Web: www.herzmed.de

The products shown are only an excerpt of our product range. Exact details and product information are available on request or at www.herzmed.de. All prices marked with an asterisk are exclusive of VAT. All brand names used remain the property of the respective brand owner. We assume no liability for errors, possible misprints and color deviations. Our general terms and conditions apply, which you can also view at www.herzmed.de.

